



2 Witney Way
Boldon Business Park
Tyne and Wear
NE35 9PE
England

Hessburg Barron Vacuum Trepine and the Barron Marking Donor Punch

The following products are now available from
Altomed Limited.

To place an order contact Customer Services on

Tel: (0) 191 519 0111 Fax: (0) 191 519 0283


The Hessburg Barron Vacuum Trepine.



This unique trephine can provide a complete or partial keratotomy for keratoplasty procedures.


Unlike conventional units, this vacuum trephine attaches to the cornea during the trephining process to hold the cornea in its natural state allowing for a perfect perpendicular cut for recipient well.

The trephine is supplied sterile and is disposable

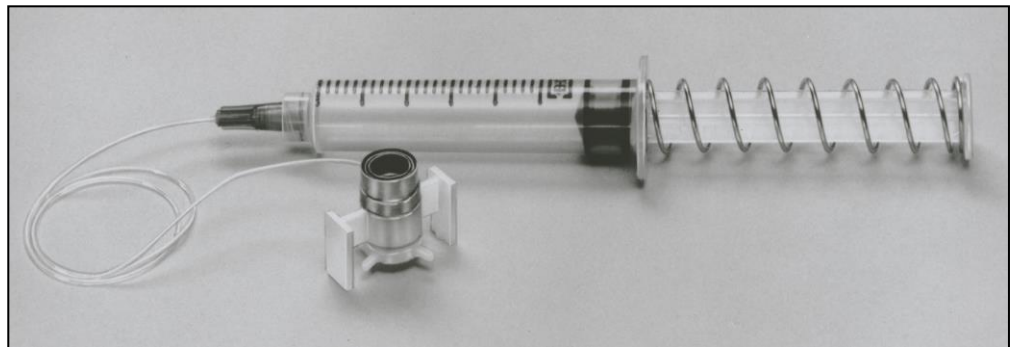

SIZE	6.0mm	6.5mm	7.0mm	7.5mm	7.75mm	8.0mm	8.5mm	9.0mm
CODE	A 7028	A 7030	A 7032	A 7034	A 7035	A 7036	A 7038	A 7040

The Barron Donor Marking Punch.

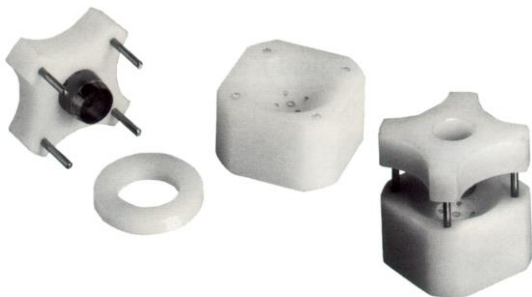

This unique disposable punch system is designed to make an exact cut of the donor tissue for corneal transplants. The ultra sharp cutting blade is guided by four posts to assure proper alignment. After the tissue is properly placed in the punch, an ink pen is used to mark the donor button in four symmetric locations. These marks will aid in the placing of primary sutures.

The system is complete with punch, blade, block and marking pen and is sterile packed ready for use. Because the system is totally disposable, you are assured a clean straight punch each time.

SIZE	6.00 mm	6.50 mm	6.75 mm	7.00 mm	7.25 mm	7.50 mm	7.75 mm	8.00 mm	8.25 mm	8.50 mm	8.75 mm	9.00 mm	9.50 mm
CODE	A 7041	A 7042	A 7043	A 7044	A 7045	A 7046	A 7047	A 7048	A 7049	A 7050	A 7051	A 7052	A 7053